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Share Info Suisse GmbH  
St.-Leonhard-Str. 35  
9000 St. Gallen  
Switzerland

### AOJ HEALTH Technology Co. Ltd.

Room 202, Building 5, Shuhe Industrial Park,  
Sanwei Community, Hangcheng Street, Bao'an District  
518126 Shenzhen, PEOPLE'S REPUBLIC OF CHINA

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UDI-DI (PPN): (01)06976313750496

**Intended use**

The insect bite healer can be used to treat insect bites or stings and is intended to relieve itching and swelling by means of heat (local hyperthermia).

**Technical data**

Dimensions:	110 x 34 x 24 mm
Model:	IBR30
Material:	ABS
Output/nominal voltage:	2W / 3V
Expected lifetime:	5 years
protection class:	IP42
Type of battery:	2 x 1.5V AAA batteries (included in delivery)
Battery life	approx. 300 treatments
Operating conditions:	10 °C – 40 °C ambient temperature, 30% – 85% relative humidity
Storage conditions:	-10 °C – 50 °C ambient temperature, 30% – 85% relative humidity

**Battery**

Type: AAA  
Chemical system: Alkaline  
Nominal voltage: 1,5 V  
Capacity: 1000 mAh  
Manufacturer: GuangDong Liwang New Energy Co., Ltd.  
Address: No. 10, LianTangJiao 2nd Road, TangXia Town, DongGuan City, GuangDong Province, China  
Home page: www.liwangbattery.com  
Recommended fire extinguishing agent: water fire extinguisher

**Safety instructions**

Please read the following points through carefully and follow the instructions before using the device.

**Danger**

- ! Choking hazard! Keep the battery out of the reach of children.
- ! If batteries have leaked, remove using chemical resistant protective gloves and protective glasses.

**Warning**

- ! Check the unit for visible damage before each use. If you detect any damage, do not use the unit.
- ! Avoid extreme heat and direct sunlight.
- ! Do not subject the device to mechanical shocks.
- ! With the exception of the cover on the battery compartment, the device must not be opened as this can lead to electric shocks. If a device is malfunctioning, please do NOT try to fix the problem yourself. Please contact the retailer in this case.
- ! Do not submerge the device in water, other fluids or flammable detergents, in order to avoid the risk of an electric shock and/or the danger of fire.
- ! The unit is not waterproof and must therefore not be submerged in water or directly exposed to any liquids.
- ! The device may only be operated by an adult.
- ! Do not use the device on sensitive areas of skin or near the eyes, mouth or mucous membranes.
- ! Do not use the device in the event of skin diseases, injuries, inflammation, redness, rashes, burns, frostbite, wounds or scars!
- ! Do not use the device in combination with ointments or creams.
- ! Consult a doctor before using the device and seek medical advice if the symptoms persist or worsen.
- ! Use of the device does not replace treatment by a doctor.
- ! If the device does not work or if pain occurs, use should be stopped immediately.
- ! All serious incidents that occur in connection with the product must be reported to the manufacturer and the competent authority of the Member State in which the user resides.


**Caution**

- ! Please use the item only for its intended purpose.
- ! The device may only be used on children aged 3 years and over.
- ! Only use the device for private use, it is not intended for commercial use.
- ! This product is not a toy! Do not allow children to play with the device.
- ! Remove all packaging material before use.

**Scope of delivery**

1x insect bite healer, 2x AAA alkaline batteries

**Device description**

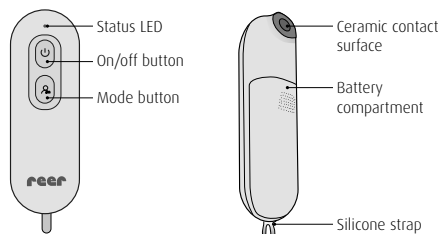**Before first use**

The batteries must be inserted into the insect bite healer for initial operation. To do this, open the battery compartment on the underside of the device. Insert the batteries and close the battery compartment again.

**Please note:** Batteries of different types or new and used batteries must not be used together. Batteries must be inserted with the correct polarity, see battery compartment.

**Functions**

If there is still a sting stuck in the skin, remove it before applying the insect bite healer.

To switch on the device, press the On/Off button on the front of the device. After a short beep, the status LED lights up green and the insect bite healer is ready for use. For best results, use the insect bite healer immediately after the bite or sting. If the time between the sting and use of the bite healer is too long, this may impair the effectiveness of the bite healer.

By pressing the mode button, you can choose between the two modes of the insect bite healer:

**Mode 1**

This mode is recommended for first-time use and for sensitive skin.

Place the insect bite healer on the skin so that the ceramic contact surface is touching the affected area. Then press the mode button **once**. The status LED will change colour to blue.

The ceramic contact surface heats up the affected area after a short heating period for 3 seconds at 47°C ±2°C (116.6°F ±3.6°F). The end of the treatment is confirmed by a short beep.

**Mode 2**

This mode is recommended for repeated use and for normal skin.

Place the insect bite healer on the skin so that the ceramic contact surface is touching the affected area. Then press the mode button **twice**. The status LED will change colour to yellow."

The ceramic contact surface heats up the affected area after a short heating period for 6 seconds at 50°C ±2°C (122°F ±3.6°F). The end of the treatment is confirmed by a short beep.

**Important!**

Depending on the sensitivity of the user and of the skin area, the treatment temperature may be perceived to be unpleasant and lead to slight reddening of the skin and, in rare cases, to skin irritation. If the heat function is perceived to be too hot, the treatment should be stopped immediately. If necessary, you can apply the insect bite healer to the affected area again after a break of 2 minutes. However, the insect bite healer should not be used more than 5 times per hour on the same area.

If the battery is too low, the status LED lights up 3 times and an acoustic signal sounds. Both batteries should then be replaced.

Switch off the insect bite healer after use with the On/Off button. The device also switches off automatically after 5 minutes without use. Before putting the device away, make sure that the ceramic contact surface has cooled down.

**Cleaning and care advice**

It is recommended that the device is cleaned before initial use and after each further use. Switch off the device and remove the batteries before cleaning. To clean the ceramic contact surface, we recommend using commercially available alcohol disinfectant wipes. Clean the housing with a soft, dry cloth only.

The device is not waterproof and so must not be immersed in water or exposed to direct moisture. Remove the batteries from the device if it is not going to be used again for a while. This product does not require special maintenance. Repairs to the device may only be carried out by a qualified professional. Store the insect bite healer in a dry, dust-free place where it is protected from dirt and direct sunlight.

**Warranty**

Please retain your receipt, the user manual and item number. Warranty claims or claims from promises of guarantee must be processed by the dealer. The technical specification, information and characteristics of the product described herein have been compiled to the best of our knowledge and belief and are correct at the time of printing. Any information provided by reer should be correct and authoritative. Nevertheless, reer does not accept liability for potential errors in this document and reserves the right to make changes to the product design and/or specifications without prior notice. Reproduction and disclosure of information contained in this copyrightprotected document, in any form or by any means - graphical, electronic or mechanical by photocopying, recording on tape or storage in a data retrieval system (even in part), require prior written consent from reer.

**Environmental protection**

At the end of its useful life, this product must not be disposed of with normal household waste. According to Section 10 (1) ElektroG, owners of old appliances are obliged to hand them in at a collection point for the recycling of electrical appliances. Please contact your local council to find out more about an appropriate disposal centre. According to Section 17 (1) and (2) ElektroG, distributors with a sales area for electrical and electronic equipment of at least 400 square metres and distributors of foodstuffs with a total sales area of at least 800 square metres are obliged to:

- a) take back the corresponding old appliance free of charge when a new appliance is supplied
- b) generally take back small waste appliances that are not larger than 25 cm in any external dimension free of charge, but limited to three waste appliances per appliance type.

Batteries and battery packs must be removed from the device before handing it over, provided that this is possible without causing damage, and separately disposed of at a collection point for batteries and battery packs.

The crossed-out wheelee bin symbol on the product, in the instructions for use or on the packaging indicates this responsibility.

**Explanation of symbols**

The following symbols are used in the operating instructions, on the packaging and on the rating plate of the appliance and accessories:

	Follow the instructions for use
	Applied part type BF
	Disposal in accordance with the Waste Electrical and Electronic Equipment Directive (WEEE)
	Manufacturer
	Date of manufacture
	Batch designation
	Protection against penetration of solid bodies with a diameter >1.0 mm; protection against dripping water with 15° inclination
	Authorised EU representative
	The CE marking certifies conformity with the essential requirements of Regulation (EC) 2017/745 for medical devices.
	Medical device
	Unique device identifier

**EU REP** Company Address

Share Info GmbH  
Am Schuzentrum 12  
41564 Kaarst  
Germany

**CH REP** Company Address

Share Info Suisse GmbH  
St.-Leonhard-Str. 35  
9000 St. Gallen  
Switzerland

**AOJ HEALTH Technology Co. Ltd.**

Room 202, Building 5, Shuhe Industrial Park, Sanwei Community, Hangcheng Street, Bao'an District 518126 Shenzhen, PEOPLE'S REPUBLIC OF CHINA

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reer GmbH