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Children-Thermometer SC 44 flex with flexible measuring point

CE 0197

Special features:

The SC 44 flex is a fast and waterproof maximum thermometer, i.e. the highest temperature is displayed after 30 seconds. The most recent measured value is always saved.

CAUTION:
A fever thermometer is not a toy.

Safety information:

- ◆ The SC 44 flex is suitable only for measuring human body temperatures.
- ◆ The thermometer contains small parts (batteries, etc) which could be swallowed by small children. Please keep the thermometer and its batteries out of the reach of small children.
- ◆ Protect the SC 44 flex from direct sunshine and temperatures above 50°C.
- ◆ Never drop the SC 44 flex, do not bend or open it.

Operating manual:

Press the I/O button to switch on. A short signal tone shows that the device is in the "switched on" condition. A segment test takes place simultaneously (checks all display segments and carries out an internal function test). After this, if the temperature measured by the temperature sensor is less than 32°C, the display shows an "L" and "C". After around one second, the "C" starts to blink. The "L" stands for low (less than 32°C). The current temperature is continuously displayed during measurement, and the "C" blinks during the entire measurement time. After the measurement is completed, the temperature value measured remains in the display until the SC 44 flex is switched off. To turn the device off after use, press the I/O button briefly. If the device is not switched off, a signal tone is sounded after around 8 minutes and the device then switches off automatically.

Fever alarm:

If the temperature register higher than 37,5°C, the fever alarm (a quick beeping sound which occurs continuously for about 8 seconds) will sound when the peak temperature has registered.

Display of last measured value:

If you keep the I/O button pressed for longer than three seconds, the last measured value will be shown in the display. As soon as you release the I/O button, the device is ready to take measurements.

Measured value deviation:

The measured temperatures will always deviate, depending on the measurement location. This is normal, and caused by deviations in temperature distribution on the human body. You could, for example, measure a rectal temperature of 36,6°C, an oral temperature of 36,1°C and an underarm value of 35,6°C.

Measurement locations:

In the mouth (oral)


Small temperature variations in the mouth between measurements are normal. These are caused by varying temperature distribution within the mouth. Try to place the thermometer in one of the two heat pockets under the tongue, to the left and right of the tongue root. The measurement sensor must have good contact with the tissue. The subject should close the mouth and breathe normally through the nose so that the measured result is not affected by breathed air.

Measuring period: approx. 30 seconds

In the anus (rectal)

This is the safest method. It is particularly suitable for infants and small children. Guide the thermometer point carefully 2 - 3 cm into the anus.

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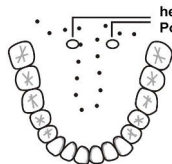
Temps de mesure: env. 30 secondes.

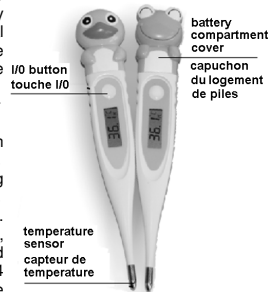
Under the arm (axillary)

From a medical viewpoint this type of measurement leads to inaccurate values.

Changing the battery:

As soon as the  symbol appears in the display, the battery has run out and requires replacement. Always replace the discharged battery with one of the same type. To do this, proceed as follows:

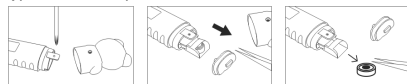Put the new battery (LR 41 or SR 41) in. Please make sure it is the right way round. You can purchase batteries from jewellers or electrical shops.

Disposal note:

Batteries must be disposed of properly. Battery disposal containers are provided in all shops selling batteries in addition to communal waste collection points. After the end of your device's service life, please dispose of the device correctly using the communal collection point for old electrical equipment.

Measurement checks

Checking the accuracy of the SC 44 flex is not necessary in most cases. The German Medical Product User regulations of 29th June 1998 stipulate a measurement check every two years for devices used in medicinal practices. If you require a measurement check (at your expense), please send the SC 44 flex to our customer service address. Please make sure it is well-packaged and that the correct postage stamps are applied.

Cleaning and disinfection:

Please only use an alcohol pad. Each pharmacy sells alcohol pads.

We recommend the use of Mucocit®-F2000 or Helipur®Hplus to disinfect the device during use in medical facilities. Both are effective against all hospital-typical pathogens and can be purchased from any dispensing chemist.

Warranty:

We give a guarantee period of 2 years from the date of purchase. We guarantee perfect design of this article.

We will rectify all material or manufacturing errors during the guarantee period. If, contrary to expectations, defects are found, please send the device back to our customer service department. Please enclose the proof of purchase and a description of the defect so that the problem can be rectified as soon as possible.

The guarantee excludes damages caused by improper handling, related to the battery or to the packaging.

Referenced norms:

The thermometer has been certified in accordance with DIRECTIVE 93/42/EU of 14th June 1993 relating to medicinal products.

The thermometer SC 44 flex complies with the requirements of European norm EN 12470-3:2000.

Specifications:

Measuring range: 32.0°C to 43.9°C
Temperature below 32.0°C: display "L" for low.
Temperature above 43.9°C: display "H" for high.
Measurement accuracy: 0.1°C between 34.0°C and 42.0°C at an ambient temperature of 18°C to 28°C
Self-test: Automatic internal check of test value 37.0°C. In case of deviation > 0.1°C, "Err" will be displayed.
Fever alarm: > 37.5°C
Storage temperature: -10°C to 50°C

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