reer

BABYWAAGE **comfort deluxe**

ArtNr: 6409

Gebrauchsanleitung Babywaage „comfort deluxe“
Instruction Manual Babyscale „comfort deluxe“
Mode d'emploi Pèse-bébé „comfort deluxe“
Istruzioni d'uso Bilancia per baby „comfort deluxe“
Gebruiksaanwijzing Babyweegschaal „comfort deluxe“
Manual del Usuario del Báscula para bebés „comfort deluxe“

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Instructions Babyscale „comfort deluxe“, Art.Nr. 6409

Congratulations on your purchase of this quality product by Reer. In order to profit from all the advantages this device has to offer, please read these instructions carefully before first use and keep them for future reference. Comply with the safety indications, so that the baby scales can reliably serve you for years.

Safety instructions / Specific function

- The baby scales are exclusively adapted for the measurement of the body weight of babies and small children up to a maximal weight of 30 kg.
- Do not overload the device. The maximal weight allowed is 30 kg. In case of overload the measuring sensor can be damaged.
- Protect the baby scales from the direct sunlight, temperatures over 50°C and the penetration of water.
- The battery included must not be charged or reactivated by other means or disassembled, thrown into the fire or shorted-out.
- Do not drop the baby scales.
- Avoid children jumping on the scales and do not open them (apart from the battery case to change the batteries).
- Avoid exposing to strong vibrations.
- Change the batteries as soon as the battery control notice appears on the display.
- Protect the scales from moisture. Use only soft cloths for cleaning.

Device description:

Technical data:

Measuring range:	0,00 to 30,00 kg
Accuracy of measurement:	0,00 - 30,00 kg + 0,02 kg (up to 5,00 kg) + 0,03 kg (from 5,01 kg to 10,00 kg) + 0,04 kg (from 10,01 kg to 15,00 kg) + 0,05 kg (from 15,01 kg to 20,00 kg) + 0,06 kg (from 20,01 kg to 25,00 kg) + 0,07 kg (from 25,01 kg to 30,00 kg)
Display:	LCD-Display
Repartition:	0,01 kg
Operational voltage:	6 Volt
Current alimentation:	4 x AA batteries (included in delivery)
Measuring unit:	kg / lb&oz
Maximal operational weight:	30 kg

HOLD FUNCTION:

As babies are always a bit agitated on the scales, the device automatically determines the value that is constant for the longest duration. So wrong values caused by touching the scale pan, or by the child kicking his/her legs, etc. can be avoided.

MEMORY FUNCTION:

After the starting up the memory function shows you the measuring value of the last measurement.

TARE FUNCTION:

With the TARE button you can implement an "additional weight function". If you want to put your baby on a pillow in the scale pan and you just want to determine the weight without the pillow, read the instructions described in paragraph 4.1. PAY ATTENTION: The sum of all the single items in the scale pan must not exceed 30 kg.

1. Pedestals	4. Music „On/Off“	7. Weighing surface
2. Speaker	5. On/off switch	
3. Battery compartment	6. LCD Display	

Instruction manual:

First starting-up:

1.1. Opening the battery compartments' cover


Open the battery compartment on the rear side of the device.

1.2. Inserting the batteries

Insert the batteries (4 pieces; type AA 1,5V) into the battery compartment. Please always pay attention to the correct polarity according to the polarity information + and - indicated in the battery compartment.

1.3. Closing the battery compartments' covers

After having inserted suitable batteries, slide the battery compartment cover back onto the battery compartment until it clicks into place.

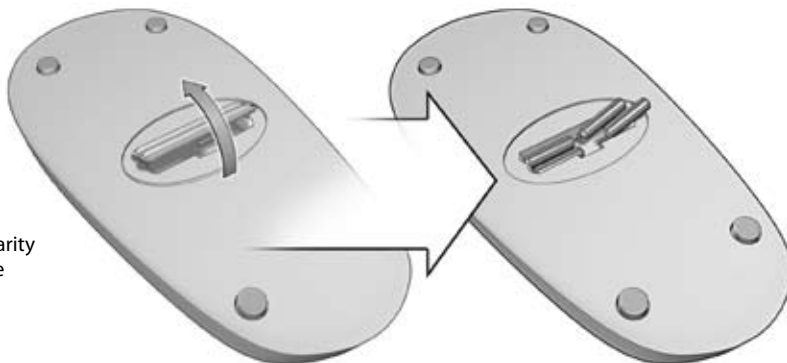2. Start-up:

The baby scales should always be placed on a smooth and dry surface. You can switch on the device pressing the button "ON/OFF"

2.1 Setting up the desired measuring unit

Set up the desired measuring unit, by pressing the "ON/OFF" button again.

The chosen measuring unit is shown on the display (kg /lb&oz).

As soon as the display changes to "0.00" (kg,lb&oz), the scales are operational and the measuring procedure can start

Memoryfunction:

At the beginning of the measurement the display indicates the measuring result of the last measurement.

As soon as the display shows "0.00", you can put the baby into the scale pan and start the measuring procedure.

3. Weighing

- 3.1 As soon as the display shows "0.00", you can put the baby into the scale pan and start the measuring procedure.

The weighing procedure is concluded, as soon as the display flashes 3 times. The measuring result is displayed on the display. As soon as you take the child from the balance, the measuring value is reset ("0.00"). This happens at the latest 30 seconds after weighing.

Automatic switch-off:

The baby scales switch off automatically after approx. 30 seconds.

4 Tare function – weighing with additional weight



If you want to carry out a measurement with additional weight (for example a towel or a pillow), the weight must be tared in order to determine the real total weight.

- 4.1 Put the additional weight (for ex. towel, pillow, etc.) on the scale frame, wait for the measuring result and press the button "On/Off" to save this value. Take the weight from the scale. The measuring result is shown on the display (image 4).
- 4.2 Continue with point 3.1

5. Start music reproduction:

To start the music reproduction, press the "Music" button

To stop it, press the "Music" button again.

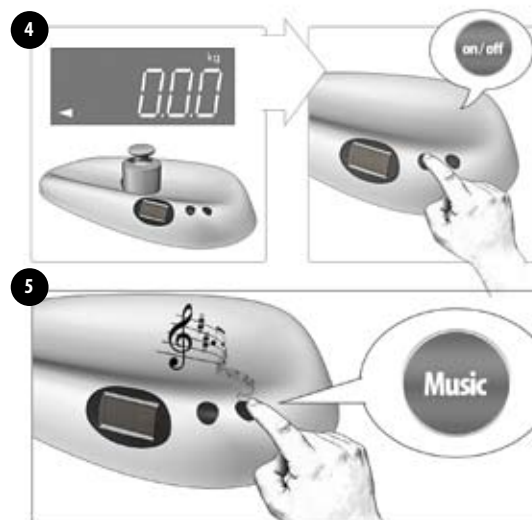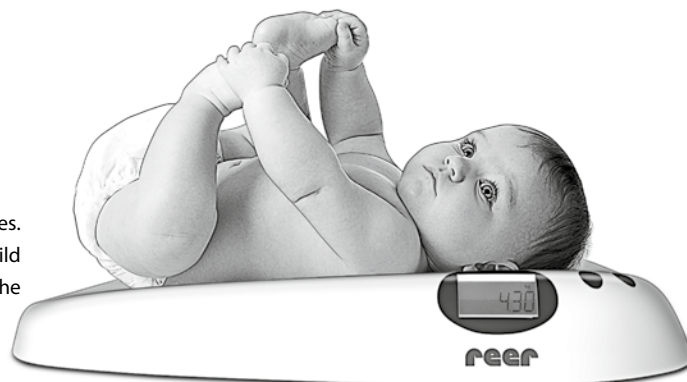Cleaning

Clean the device with a soft and dry cloth.

Never use petrol, thinners or similar corrosive solvents.

Waste disposal indications

Batteries must be correctly disposed. For this purpose the battery retailer and the municipal collection points provide suitable containers for battery disposal.

After the conclusion of the service life, please dispose the device correctly at the municipal collection point for used electrical devices

Guaranty:

Please keep a record of the purchase receipt and article number. Any warranty or guarantee based claims are to be directed solely to the retailer.

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All rights reserved, technology and design subject to change without notice.
Reer GmbH, Riedwiesenstraße 7, D-71229 Leonberg
reer GmbH, Postfach 1639, 71206 Leonberg, www.reer.de